
LAFE

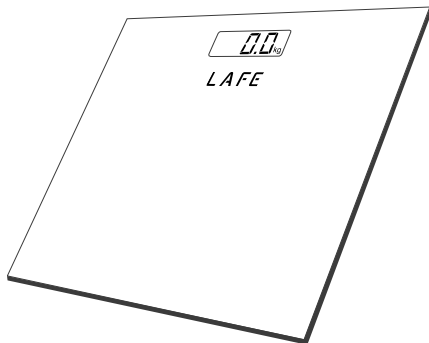FitScale

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EN Electronic bathroom scale

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Zakres pomiaru: od 3 – 150 kg

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WARNING

The device is intended for domestic use only.

The device is NOT intended for commercial/professional use.

The device is NOT intended to be used outside.

Note! Prior to use, you should familiarize yourself with this instruction manual in order to avoid accidents and to ensure proper use of the appliance. Manual should be kept and stored so that it is always easy to reach.

SAFETY INSTRUCTIONS FOR USE

1. Before first use, remove all packing elements.

Note! Housing with metal elements may include barely-visible protective film stretched over those parts, which needs to be removed.

2. The appliance should not be used by children or persons with reduced physical, sensory and mental ability without attendance of qualified or experienced persons and should always be used in accordance with the instruction manual.

WARNING: This appliance can be used by children over the age of 8 and persons with limited physical, sensory or mental capabilities, or persons without experience and knowledge of the appliance, ONLY if used under the supervision of a person responsible for their safety or if they have been given guidance on safe use of the appliance and are aware of the hazards associated with its use. Children should not play with the appliance. Cleaning and maintenance operations should not be performed by children unless they are over the age of 8 and these operations are performed under supervision.

3. NEVER PLACE ANY METAL ITEMS inside the appliance.
4. Use the device on dry and stable surface only.
5. Never jump on the device.
6. Do not use the device with your legs wet or when there is water on the device as there is a risk of slipping.
7. Remove batteries if the device is not used for a long period of time. Otherwise the batteries may leak, which may damage the device permanently.
8. The manufacturer shall not be liable for any damage caused by improper cleaning of the appliance.

Warnings

1. Never use the appliance if it is damaged or works abnormally.
2. Do not use the appliance if it has been dropped from a height and shows visible signs of damage.
3. Any repair, removal or replacement of any parts should always be performed at a specialized facility.
4. If you think the appliance might be broken, NEVER try to repair it yourself.
5. Never wash the device in running water or in any manner, which allows water to drip onto the device.
6. The appliance should not be exposed to extreme temperatures, direct sunlight, excessive humidity or excessive dust conditions.
7. Do not use the appliance near heaters, open flames and other sources of heat.
8. The appliance is intended solely for weighing people.

Basic information

1. 4 precision sensors of pressure
2. LCD Display
3. Unit switch /
4. Battery compartment

Remarks concerning use

1. The scale must be used on a hard, flat surface such as parquet, tiles, linoleum, etc. We do not recommend using the scale on carpets and soft flooring.
2. Prior to first use, remove insulation from the battery (Figure 2).
3. Prior to first use, set the appropriate units by pressing the button underneath the platform (Fig. 3) (Fig. 4).
4. When weighing, stand on the scale in a calm and stable fashion.
5. The weighing results may be distorted after intense exercise, diet or in cases of extreme dehydration.
6. We advise you to check your weight at the same time of day, preferably in the morning after waking up and before eating, exercise or other physical activity.
7. Keep your scale in a horizontal position – it is safer for the sensors.

8. When the scale is not in use, do not place items on it which may unnecessarily burden the sensors located in the feet of the scale.
9. The scale operates in the range from 3 kg to 150 kg.
10. If the display shows „Err“ (Fig. 5) – it means that the burden is too big.
11. If the display shows „Lo“ (Fig. 6), it means that the battery is low and needs to be replaced.

Operation

Power supply – Battery:

1. If the scale is new and unused, remove the safety seal first (Fig. 2).
2. If the batteries are weak or completely exhausted, replace them in accordance with the following instructions:

Note: If the batteries are running low, this may affect the correctness of the measurement.

- a) Open the cover of the battery compartment (Fig. 6) (Fig. 7).
- b) Use a thin pointed object to pry the battery out and remove it.
- c) Insert new battery in accordance with polarity indicators – plus up, and press lightly.
- d) Finally, replace the battery compartment cover.

It is recommended to solely use the disposable batteries to supply the device.

Important safety information

Equipment supplied with batteries

The following safety measures concern products supplied with rechargeable or disposable batteries.

Inappropriate use of batteries may result in electrolyte leaks or explosion. Released electrolyte causes corrosion and may be toxic. It may also cause skin and eye burns; it is also toxic when swallowed.

To reduce the risk of injury:

- Batteries should be kept out of children's reach.
- Batteries should not be heated, open, pierced, destroyed or put into fire; they should not be inserted into the device backwards. Pay special attention to the pole marking [+] and [-].
- Do not insert old and new batteries at the same time, or batteries of different types (for example carbon and zinc or alkaline batteries).
- Do not allow the contact of the metal ending of the batteries with metal objects. These elements may heat up and cause burns.
- If the batteries are used up, or if the device is to be stored for a longer period of time, remove the batteries from it.
- Old or used batteries should be removed from the device and utilised or recycled in accordance with the domestic regulations on waste disposal.
- In case of leaks from the battery, remove all batteries avoiding contact of the skin or clothes with the electrolyte. In case of contact with the electrolyte rinse the skin or clothes with water. Before installation of new batteries clean the battery compartment with damp paper towel.

Only disposable batteries

- Caution! Replacement of batteries into inadequate type may result in explosion. Use batteries of the same type and size (alkaline, zinc-carbon or zinc and chlorine).
- Do not recharge the traditional batteries

Chemical composition of the applied battery:

Lithium (Li-ion)

Mode of removal, replacement:

Follow the battery replacement procedure.

Turning the scale on and measurements (Fig. 8)

1. Place the scale on a hard, flat surface.
2. Press the platform vigorously with one foot (you can also stand on the platform and then get off).
3. The display will show numbers – which means that the scale is turned on.
4. If the scale indicates zero, it is ready to weigh (be sure to check if the units in which you want measurements to be displayed are set correctly).
5. Stand on the scale in a stable, upright position.
6. After about 5 seconds, the scale will begin to measure.
7. Once the measurement is completed, the result will begin to blink.
8. If you want to then make a second measurement you need to zero the scale by pressing vigorously on the platform – the display should indicate zero.
9. After weighing, you can set the scale aside. Then after about 5 seconds of rest, the scale turns off automatically.
10. If the scale displays random values or doesn't turn off automatically it must be „restarted“. To do this, remove the battery for about 15 seconds and put it back in. The scale should begin to work again properly.

Cleaning and maintenance

Note! Always remove the batteries before cleaning.

The appliance may be cleaned with a damp cloth.

Do not use harsh cleaning chemicals that can damage the surface.

Technical specification

Model: LAFWAG44590 / LAFWAG44591

Power supply: Battery 3V 1 x CR2032

Size / weight of the appliance: 30 x 30 x 2 cm / net weight: 1.49 kg

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The symbol of a crossed-out wheeled bin placed on electronic or electrical equipment, its packaging or accompanying documents means that the product may not be thrown out together with other waste. Used equipment may contain substances with toxic and carcinogenic properties, hazardous to human health and life, and poisoning the soil and groundwater. It is the user's responsibility to hand over the used equipment to a designated collection point for its proper processing. For more information on recycling of electronic and electrical equipment, please contact your local authorities, waste disposal services and the place where you purchased this product.

The household plays an important role in contributing to the reuse and recovery of secondary raw materials, including recycling, waste equipment. At this stage, attitudes are formed that affect the preservation of the common good, which is a clean natural environment.

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